

Supermouth® quick start guide

Revolutionizing Oral Care


1



attach the brush head

- remove brush head from the packet (if not already attached)
- attach by sliding head onto the metal rod of the SmartBrush handle
- store the additional brush head in the packet; replace in 6 weeks

2

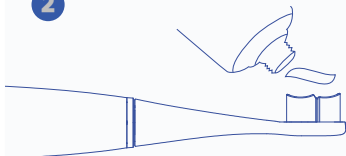


apply toothpaste

- always rinse the brush head
- apply a pea-size amount of toothpaste to the bristles
- put the brush head into your mouth before turning on the brush (to avoid spraying paste around your bathroom!)

for best results, use with SuperMouth toothpaste

3

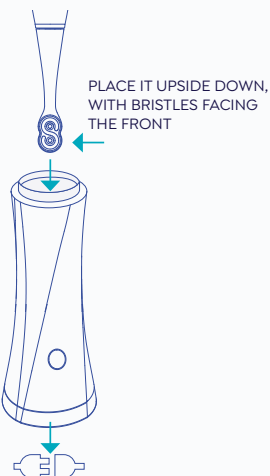

using your SmartHub to store, charge, and maintain your SmartBrush

the SmartHub is your all-in-one solution to store, charge, dry, and sanitize your SmartBrush


1. remove the clear rinse cap on the SmartHub
2. insert the SmartBrush by placing it upside down with bristles facing front
NOTE: Always remove the brush head protective cap
3. connect the USB cable to the port on the back of the SmartHub and plug into a wall outlet
4. snap the clear rinse cap back into place

light indicators

1. SmartHub: blue = UV sanitizing (5 mins), red = charging/drying.
2. SmartBrush: pulsing red = charging, solid red = fully charged

The SmartMag adapter is an alternate charging option and should only be used when the SmartHub is not available. Always use the SmartHub to store and charge your SmartBrush whenever possible.

4

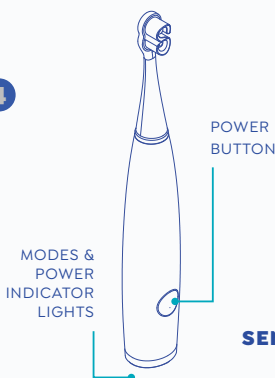

select your preferred mode

- to turn on/off your SmartBrush, press the power button
- the ULTIM8 Sonic SmartBrush has two modes designed to effectively clean your teeth
- to switch modes, turn the toothbrush on and press the button within 10 seconds; the sequence for mode switching is: first mode > second mode > off > first mode


SUPER

operates at full power to deliver a vigorous and thorough cleaning. in this mode, the bottom of the handle turns white

SENSITIVE

reduced sonic vibrations create a gentler cleaning experience for those with sensitive teeth and gums. in this mode, the bottom of the handle turns blue

5

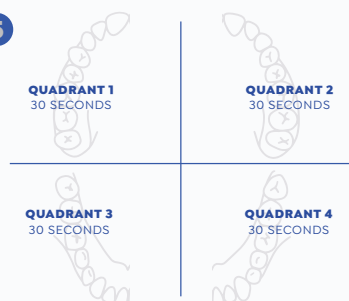

start brushing using the smart timer & pressure sensor

start with the outside surfaces of the teeth. guide the brush head slowly from tooth to tooth, holding the brush head in place for a few seconds against each tooth before moving on to the next one. following along with the shape of each tooth and the curve of the gums is easy with the brush head designed to hug every curve.

about the smart timer

every 30 seconds while brushing, you'll feel a short pause of the vibration, signaling it's time to switch to a different section of your mouth. this ensures each of the four areas receives equal attention. after two minutes, the toothbrush automatically turns off.