

Supermouth®



ultim8™

Travel SmartCase™

user manual

welcome

you deserve the ultimate
brushing experience

table of contents

introduction	1
what's in the box	2
smart features	3
getting your Ultim8 Travel SmartCase™ started	5
Ultim8 Travel SmartCase™ precautions	7
warnings	9

ultim8™

the perfect integration of form and function
for the ultimate brushing experience




The Ultim8 Travel SmartCase™ is the ultimate travel companion for your Ultim8 Sonic SmartBrush™. Designed for discerning travelers, it offers an elegant and convenient solution for carrying and protecting your toothbrush, SuperMouth travel floss, and travel toothpaste.

The case features a state-of-the-art UV sanitizer that eliminates 99.9% of germs on your toothbrush, ensuring the highest standards of hygiene. It also includes a sleek charging cable and a rechargeable battery, allowing for effortless recharging and on-the-go sanitizing, making it the perfect blend of luxury and practicality for your travels.

what's in the box



Ultim8 Travel SmartCase™

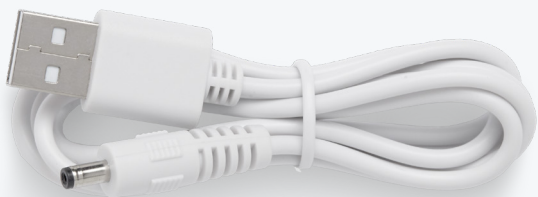

USB cord

Ultim8 Travel SmartCase™

with *8 smart features* for
unparalleled hygiene and
convenience

- 1 stylish & ergonomic design**
is compact but still holds essential
brushing items for easy travel
- 2 TSA approved lithium battery**
permits carry-on and checked
luggage for hassle-free travel
- 3 hygienic storage and recharging** of
your Ultim8 Sonic SmartBrush
- 4 UV sanitizer** effectively eliminates
99.9% of germs in 3.5 minutes
- 5 efficient charging** for your
SmartBrush using SmartMag
technology and a built-in battery that
allows the UV sanitizer to run for up to
3 weeks on a single charge
- 6 storage** for SuperMouth travel-size
toothpaste and floss (sold separately)
- 7 quick-access latch** simplifies interior
cleaning and allows for efficient
water drainage
- 8 integrated cable management**
universal charging cable compatible
with all Ultim8 products

getting started with your **Travel SmartCase**

usage and charging instructions

- 1.** prepare SmartBrush by rinsing the brush head thoroughly; shake/wipe off excess water from brush head and handle.
- 2.** remove brush head from your SmartBrush.
- 3.** place SmartBrush handle in designated slot ensuring the front faces forward.
- 4.** place brush head in SmartCase well; ensure bristles face the brush handle (see picture).
- 5.** place SuperMouth travel floss and toothpaste in designated spaces (see picture); these items are sold separately at supermouth.com.
- 6.** plug charging cable into a USB wall charger (not included); insert other end into SmartCase.
- 7.** close SmartCase lid firmly to activate the UV sanitizer which operates for 3.5 minutes and automatically shuts off when complete.
- 8.** during sanitization, a solid blue light indicates the process is in progress. Once sanitization is complete, the SmartCase light will turn solid red if plugged in or turn off if not plugged in.

charging Indicators:

a yellow light or failure to turn on indicates the need to charge the SmartCase.

a red light indicates that the SmartCase is charging.

a green light indicates that the SmartCase is fully charged.

NOTE: SmartCase must be plugged in to activate charge; will not charge if battery is fully charged

SmartCase UV sanitizer operates up to approximately 3 weeks @ 3.5 min twice a day on rechargeable lithium ion battery even when the case isn't plugged in.

WARNING: Do not touch or eye the UV light of the sanitizer when in use, UV light can be harmful to the human eye and skin. The UV light bulb of the sanitizer is hot during and immediately after the sanitizing cycle. Do not touch the UV light bulb when it is hot.



cleaning SmartCase™

1. wipe inside with microfiber cloth once per week when using.
2. lift-up water drip plate under the brush head and wipe with microfiber cloth to remove any water and/or residue.
3. make sure drain hole under the water drip plate is open and free of debris so water can easily drain.
4. replace water drip plate.
5. make sure inside is thoroughly dry before storing the SmartCase to prevent formation of mold.
6. store unplugged in a cool and dry place away from direct sunlight.



function description

- automatic recognition when the handle is placed, removed and the battery is fully charged, energy saving design is adopted
- wireless charging function
- USB interface, simple power supply, easy to use
- the charging base can be magnetically fitted with the matching toothbrush

precautions

- if the product does not work for the first time, the electrical outlet may be in poor contact, please replug into electrical outlet and try again
- stop using product if anything unusual is found
- do not place the product in a place where it is easy to fall

technical parameters

Rated working voltage: DC 5V

Rated working power: $\leq 3.0W$

Rated input voltage: DC5V

Rated input power: $\leq 3.0W$

Static power: $\leq 0.5W$

Battery type: Lithium ion battery

Battery specifications: 0.925Wh

Battery life: 3 weeks @ 3.5 min twice a day

Working environment temperature: 41 – 104°F
(5–40°C)

Waterproof: IPX2

Storage temperature: Low temperature: -22°F
(-30°C), 4 hours; High temperature: 140°F (60°C),
humidity 90%, 4 hours

Adaptive product description: for SG-2739A

Charging time: ≤ 6 hours

warnings

Always wipe off all excess water from your toothbrush prior to sterilization. This product is not intended for use by persons (including children) with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the product by a person responsible for their safety.

Cleaning shall not be performed by children. Do not clean the product in the dishwasher.

Do not use attachments other than those recommended by SuperMouth.

This product contains no user-serviceable parts. The product is not intended to be repaired by service personnel in case of failure or component defect (unit can be thrown away).

Do not open and disassemble the device. Discontinue use if the product appears damaged in any way. This product complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer prior to use.



This appliance and its accessories comply with part 15 of the Federal Communication Commission regulations. Operation is subject to the following two conditions; (1) this device may not cause harmful interference, and (2) this device must accept any interference received; including interference that may cause undesired operation.

The End user must follow the specific operating instructions for satisfying RF exposure compliance.

danger

In order to reduce the risk of electrocution:

1. Do not place or store the travel case or charger where it can fall or be pulled into a tub or sink.
2. Do not place the travel case or charger in water or other liquid.
3. Do not reach for a travel case or charger that has fallen into water. Unplug immediately.
4. Do not use it while bathing.

warning


To reduce the risk of burns, electrocution, fire, or injury:

1. This product uses a built-in Lithium battery. To avoid overheating, fire or explosion, keep away from extreme heat. Before discarding, remove the battery and safely dispose of it.
2. Close supervision is necessary when this product is used by, on, or near children or invalids.
3. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by SuperMouth.
4. Keep the cord away from heated surfaces.
5. Never drop the appliance or insert any object into any opening of the appliance.
6. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged or if it has been dropped into water.

important

1. Always shake or wipe off excess water from your toothbrush and bristles before storing it in the case to help prevent mold.

2. Do NOT use the device while charging.
3. Do NOT expose your eyes or skin directly to UV-C light.
4. Do NOT drop the device or immerse it into water.
5. Always shake off all excess water from your toothbrush prior to sterilization.
6. Do not clean the product in the dishwasher.
7. Do not touch the power plug with wet hands. This can cause an electric shock.
8. When unplugging, always hold the power plug instead of the cord.
9. Do not modify or repair the product. This may cause fire, electric shock, or injury. Do not disassemble the product except when disposing of the battery.
10. This product is not intended for use by children under age three.

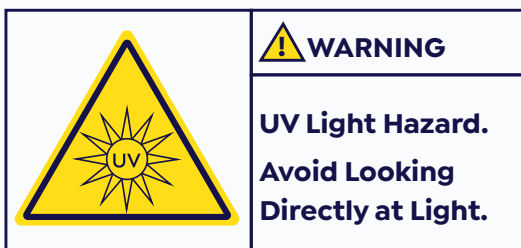

Ultraviolet radiation. Follow instructions.

DANGER-Ultraviolet radiation. Avoid overexposure.

As with natural sunlight, overexposure can cause eye and skin injury and allergic reactions. Repeated exposure may cause premature aging of the skin and skin cancer. You **MUST** avoid exposure to direct or reflected germicidal ultraviolet rays.

CAUTION

If you are undergoing treatment for any oral care condition, consult your dental professional prior to use. SAVE THESE INSTRUCTIONS.

ENVIRONMENTAL NOTICE


This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Please, dispose of the product at an approved recycling or disposing location in accordance with local regulations.

RECYCLING



This product shall not be disposed of with normal household waste.

Please take your product to an official collection point or to a service center to have a professional remove the rechargeable battery.

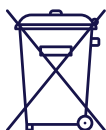

This product contains a built-in rechargeable battery which shall not be disposed of with normal household waste.

The UV light bulb of the sanitizer contains mercury. Do not throw the UV light bulb away with the normal household waste at the end of its life. Hand it in at an official collection point for recycling.

Follow your country's rules for the separate collection of electrical and electronic products and rechargeable batteries. Correct disposal helps prevent negative consequences for the environment and human health.

DISPOSAL



This product contains a rechargeable battery lithium-ion battery which must be disposed of properly.

Contact your local town or city officials for battery disposal information. You can also call 800 822 8837 or visit www.call2recycle.org for battery drop-off locations.

EXCLUSION FROM WARRANTY

The following items are not covered under the one year warranty,

- Damage caused by use of unauthorized replacement parts or unauthorized products
- Damage caused by misuse, abuse, neglect, alterations, or unauthorized repairs
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading

Changes or modifications not expressly approved by SuperMouth could void the user's authority to operate the equipment.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television

reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna
- Increase the separation between the equipment and receiver
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected
- Consult the dealer or an experienced radio/TV technician for help

.....

CANADA:

This device contains license-exempt transmitter(s)/ receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

- This device may not cause interference.
- This device must accept any interference, including interference that may cause undesired operation of the device

L'émetteur/récepteur exempt de licence contenu dans le présent appareil est conforme aux CNR d'Innovation, Sciences et Développement économique Canada applicables aux appareils radio exempts de licence. L'exploitation est autorisée aux deux conditions suivantes:

- L'appareil ne doit pas produire de brouillage;
- L'appareil doit accepter tout brouillage radioélectrique subi, même si le brouillage est susceptible d'en compromettre le fonctionnement

Supermouth®

supermouth.com

27000062-00